

TO WIN UNIVERSAL CONFIDENCE WITH KDL QUALITY AND REPUTATION
创KDL品牌 取信亿万用户



上海康德莱集团 Shanghai Kindly Group

南昌康德莱医疗科技有限公司

Nanchang Kindly Medical Technology Co.,Ltd

地址:江西省南昌市高新区艾溪湖北路新城吾悦写字楼A座7楼712室

(Add):F/7,Tower A ,No.77 North Aixihu Road,Nanchang High-Tec Industry Development Zone,
Nanchang,Jiangxi 330096.China

邮编(Postcode): 330000

电话(Tel):86-791-86861216

传真(Fax):86-791-86861216

www.kdInc.com

AMOKAN
爱默康

Veterinary Health Catalogue 动物医疗



Veterinary Hypodermic Needles 兽用注射针

Veterinarians use disposable needles to inject animals. But that always can't meet the requirement of the connecting strength and the rigid due to the particularity of the animals. Because the needles may stay in animals, and the meat with needle will hurt people. So we must use the special veterinary hypodermic needle for animal injecting.

兽医在使用一次性针头给动物注射时，由于动物的特殊性，普通针头往往不能满足连接强度和刚性的要求。导致针可能留在动物体内，而带针的肉制品会伤害到人类。因此，我们必须使用兽用注射针进行注射。




Product features 产品特点



SUS 304 cannula connected to needle hub by aluminum rivet reinforcement treatment can effectively prevent falling off. 304 stainless steel through aluminum rivet connection to the needle hub, high connection strength can avoid falling off.



Needle of pen type is much safer for veterinary injecting, and more convenient to carry. 笔式护套设计方便运输和携带。

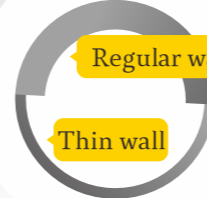


Regular wall needles are less likely to bend. 正常壁设计，防止折断。



Color-coded Poly hub for easy gauge identification. 针座带有色标，便于识别规格。

Veterinary 兽用注射针 Hypodermic Needles

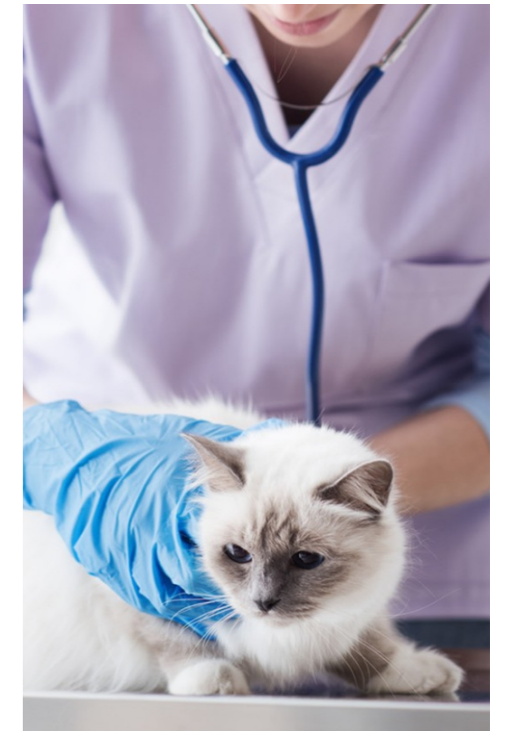

Veterinary Hypodermic Needles 规格型号表

Imperial Size 英制规格	Metric Size 公制规格	Length 长度
30G	0.3	13 ≤ L ≤ 50
27G	0.4	
26G	0.45	
25G	0.5	
24G	0.55	
23G	0.6	
22G	0.7	
21G	0.8	
20G	0.9	
19G	1.1	
18G	1.2	

Veterinary Hypodermic Needles (Aluminum Hub) 兽用注射针(铝座)

Veterinarians use disposable needles to inject animals. But that always can't meet the requirement of the connecting strength and the rigid due to the particularity of the animals. Because the needles may stay in animals, and the meat with needle will hurt people. So we must use the special veterinary hypodermic needle for animal injecting.


兽医在使用一次性针头给动物注射时，由于动物的特殊性，普通针头往往不能满足连接强度和刚性的要求。导致针可能留在动物体内，而带针的肉制品会伤害到人类。因此，我们必须使用兽用注射针进行注射。



Veterinary Hypodermic Needles 兽用注射针



Aluminum hub for large animal applications requiring a durable needle.
采用铝制针座可以在大型动物注射时更加耐用。

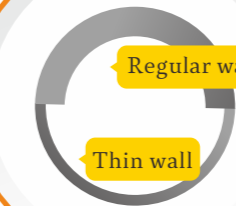

Regular wall needles is less likely to bend.
正常壁设计，防止折断。



Tri-bevel point, siliconized for smooth penetration.
三斜面针尖，硅化后穿刺更锋利。




Needle of pen type is much safer for veterinary injecting, and more convenient to carry.
笔式护套设计方便运输和携带。

Veterinary Hypodermic Needles
规格型号表

Imperial Size 英制规格	Metric Size 公制规格	Length 长度
23G	0.6	13 ≤ L ≤ 50
22G	0.7	
21G	0.8	
20G	0.9	
19G	1.1	
18G	1.2	
16G	1.6	
14G	2.1	

Veterinary Hypodermic Needle (Cooper Hub) 兽用注射针 (铜座)

Veterinary Hypodermic Needles is Ultra-sharp, tri-bevelled and Sterile needle, which is best suitable for veterinary use.
兽用注射针非常锋利, 三斜面无菌注射针适合兽医使用。

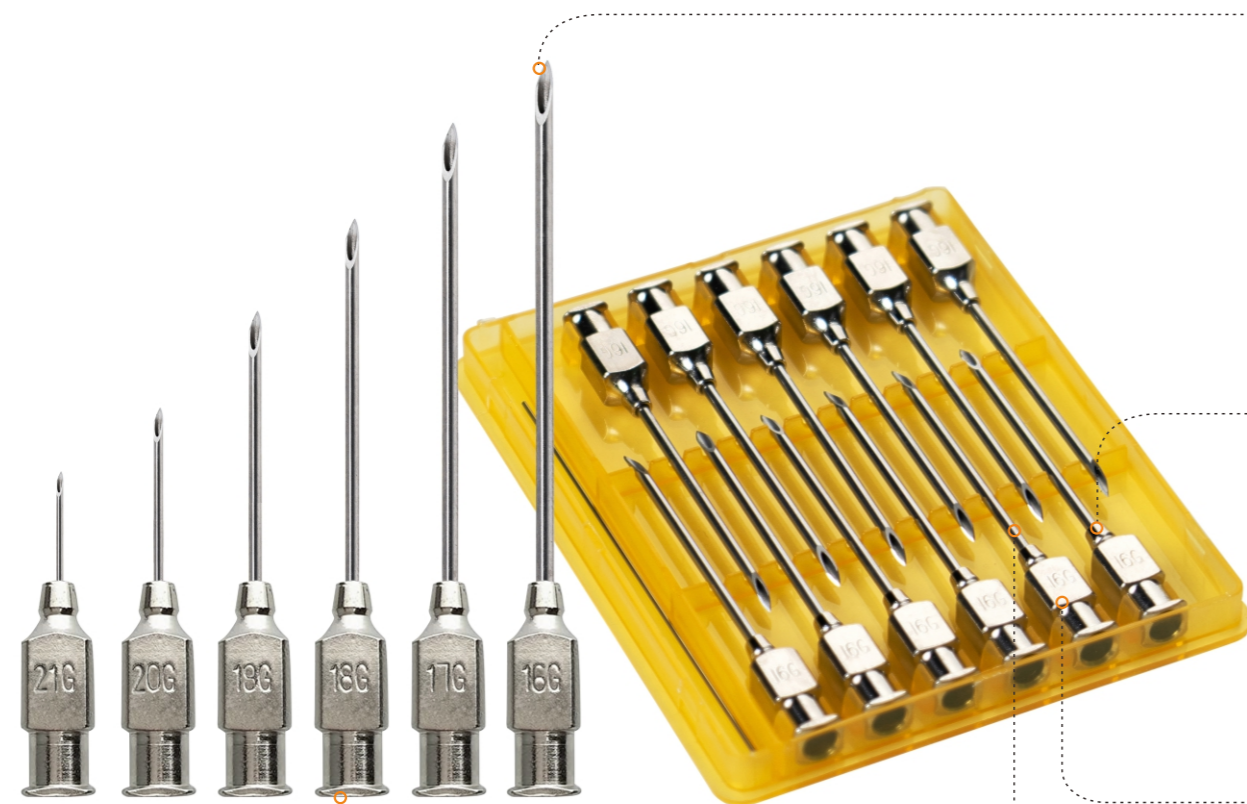

Luer-Lock

- Luer-Lock can be available in square and round Hub and HUB is made from nickle plated brass.
镀镍黄铜材质的鲁尔接口针座。
- Leak proof joint between Hub and cannula prevents cannula come out of Hub during injection.
针座和针管之间密闭设计, 防止注射过程中针管脱落。



Cannula made of stainless steel surgical grade steel, triple bevel sharp point grinding for easy penetration.
三斜面不锈钢针尖设计, 穿刺更锋利。



The SUS 304 cannula is pressed directly with the nickel plated copper material base.
304不锈钢针管与铜座压合, 连接强度高。



Stamping marks on the hub is for easy identification of needle size.
规格标记在铜座表面, 清晰可辨。



Thick walled cannula prevents needle point bending during repeated use.
厚壁针管, 防止折断。

Veterinary 兽用注射针 (铜座) Hypodermic Needles

Veterinary Hypodermic Needles 规格型号表

Imperial Size 英制规格	Metric Size 公制规格	Length 长度
20G	0.9	13 ≤ L ≤ 50
19G	1.1	
18G	1.2	
16G	1.6	
15G	1.8	
14G	2.1	



Veterinary IV Catheter 兽用留置针

- PU tubing material. PU可软化导管。
- Small and fixed wings for increased patient comfort. 小翼片更舒适。
- Thin-walled catheter with a large inner diameter assures good flow. 薄壁侧孔针管, 回血更快。
- Available in the popular size 26G for small species, including reptiles and birds. 26G更加适合于小动物。



Gauge 规格	Catheter Length (mm) 导管长度	Outer diameter (mm) 外径	Flow Rate (ml/min) 流速	Color Code 色标	Shelf Box / Shipping Case 装箱数量
26G	15	0.6	12	Purple ●	100/2000
24G	19	0.7	22	Yellow ●	
22G	25	0.9	36	Blue ●	
20G	32	1	58	Pink ●	
18G	32	1.2	90	Green ●	
17G	45	1.5	120	White ○	
16G	45	1.8	172	Grey ●	
14G	45	2	208	Orange ●	

Product features



Radiopaque lines for
X-ray visibility

X射线下可视导管



The mini wings
for patient comfort

小翼片贴合皮肤更舒适



Colour-coded hub for easy
size recognition

导管座有色变便于
识别规格




Flash-back window
可视回血窗

可视回血窗



Isolation of blood
隔离塞防止血液溢出

隔离塞防止血液溢出

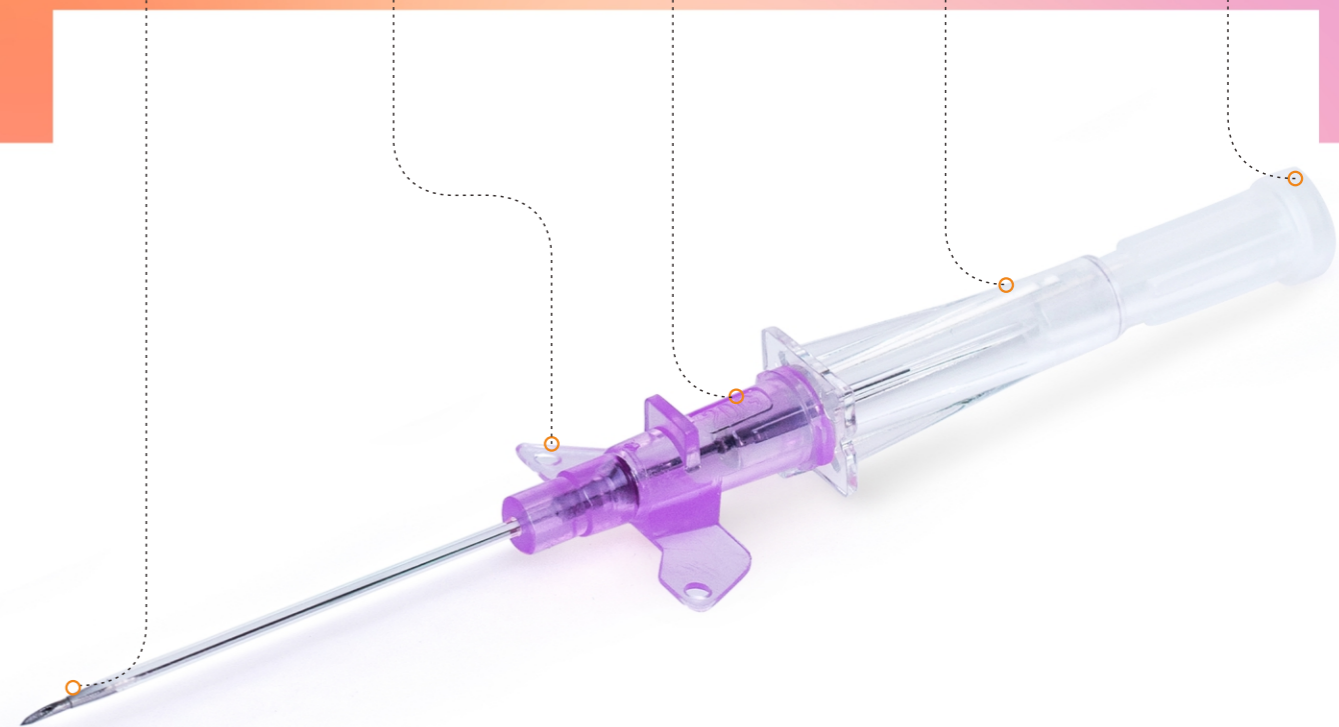

Back Cut Bevel
Extremely sharp
and back-cut for
smooth Insertion
反倒角针尖穿刺锋利



PU Biomaterial
Catheter
导管可在血液中软化,
保护血管

Veterinary Urethral Catheter 兽用导尿管

- Catheter is made of Low-density polyethylene (LDPE). 导管由低密度聚乙烯材料 (LDPE)制成。
- The catheter is rounded ends and 2 atraumatic lateral holes. 导管末端圆头并带有两个侧孔。
- Our Catheters has been expanded to include new features, materials and sizes to cover the different needs in the feline patient. 导管可以应用于各种大小动物。



Rounded ends to avoid tissue damage
圆头设计, 防止损伤。



PE material, Soft & Smooth
PE材料, 柔软舒适。



metal rivet reinforcement treatment can effective prevent falling off
导管和管座通过铆钉连接, 牢固不脱落。



Period of placement 1-3 days
留置时间为1-3天

Veterinary Urethral Catheter 导尿管规格

Type Size 种类	Metric Size 公制规格	Imperial Size 英制规格	package 装量
Open end 开放式	1.0×130mm	3Fr×5.2in	12Pcs(支)/Box(盒)
	1.3×130mm	4Fr×5.2in	
	1.0×100mm	3Fr×3.9in	
Side holes 侧孔	1.0×130mm	3Fr×5.2in	
	1.3×130mm	4Fr×5.2in	
	1.3×130mm	3Fr×6.3in	
	1.0×130mm	3Fr×5.2in	
	1.3×130mm	4Fr×5.2in	
	1.0×160mm	3Fr×6.3in	

Veterinary Easy Slide Catheter 兽用超顺滑导管

- Catheter is made from DEHP-free medical grade of PTFE. 导管由不含DEHP的PTFE材料制成。
- The catheter is soft and smooth, which ensures easy catheterization. 导管柔软顺滑, 方便使用。
- The suturing plate is soft and with suturing holes, fitting with the syringe by luer lock. 圆形缝合板柔软, 方便固定, 鲁尔圆锥接口可与注射器直连。



Rounded ends to avoid tissue damage
圆头设计, 防止损伤。



Suitable for medium-long placement (3 days)
适用于中长时间留置 (3天)

Lure-lock fitting 鲁尔接口

Handle 手柄

Suture holes 缝线
secure catheter to the prepuce
缝合孔方便将导管固定在包皮上

metal rivet reinforcement treatment can effective prevent falling off
金属铆钉连接, 防止脱落

PTFE material, Soft & Smooth
PTFE材料, 柔软舒适

Veterinary Easy Slide Catheter 超顺滑导尿管规格

Type Size 种类	Metric Size 公制规格	Imperial Size 英制规格	package 装量
Side holes 侧孔	1.2×110mm	3.5Fr×4 1/3 in	5Pcs(支)/Box(盒)
	1.2×140mm	3.5Fr×5 1/2 in	
Open end 开放式	1.2×110mm	3.5Fr×4 1/3 in	
	1.2×140mm	3.5Fr×5 1/2 in	
	1.2×180mm	3.5Fr×7 in	